



ELI  
**EZER**

# EDC-2600

Digital Chart LCD

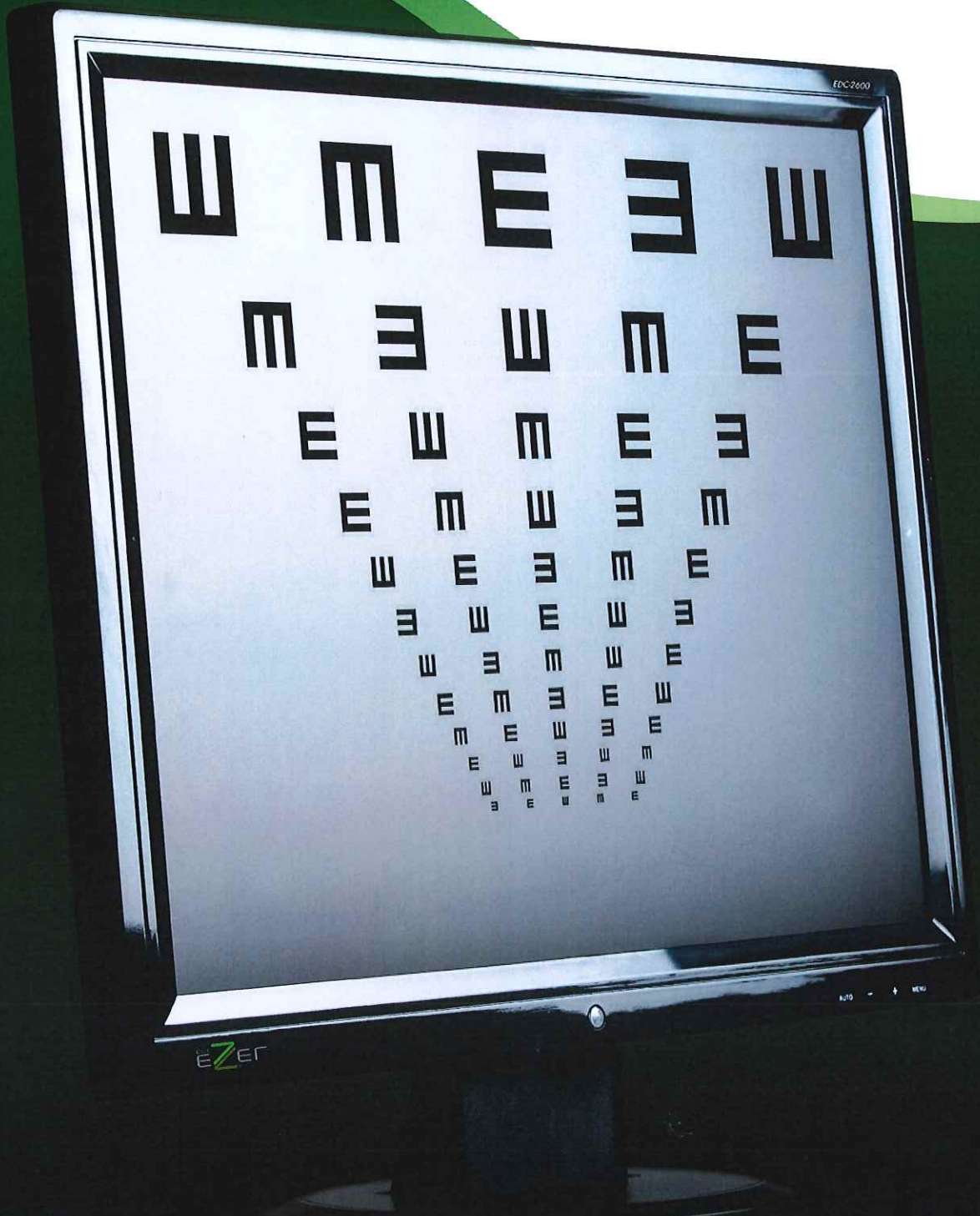


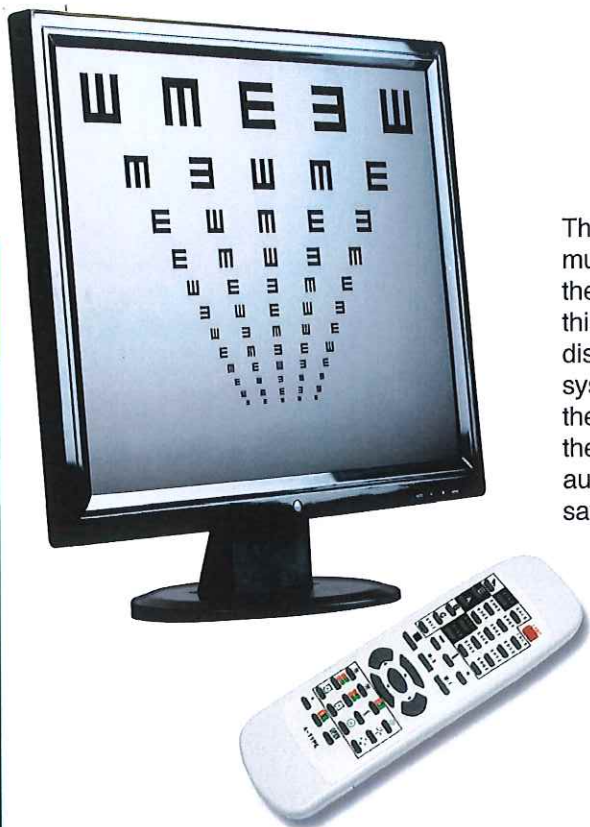

The EDC-2600 is a 19" digital chart LCD screen that is a must have for the office. There will be no need to move from the refractor to the projector during the refraction. In addition, this projector will automatically calibrate to the working distance in your office, thereby eliminating the need for mirror systems or technicians. Patient education is very important, therefore "The Eye" diagram is essential when explaining to the patients the findings discovered during the exam. The auto-saving mode kicks in after 10 minutes of inactivity saving you electricity.

### The EDC-2600 19" Panel

A high quality and outstanding clarity chart that will simplify the examination process.

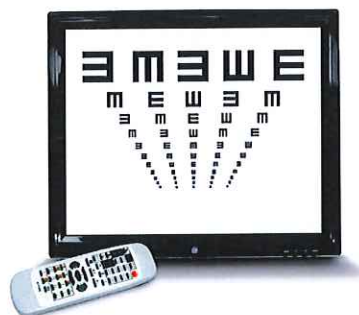


### Vertical & Horizontal Line Isolation

This feature will allow you to show the patients exactly which line or letter you want them concentrate on, thereby providing you with more accurate results.




### "The Eye" Diagram

Whenever there is a need to explain or demonstrate something to a patient, you can easily bring up "The Eye" diagram. This will allow the patient for a better understanding of the findings.

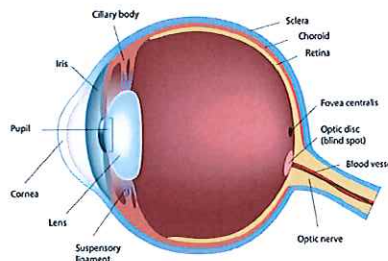

### Multiple Installation Options

The EDC-2600 can be wall mounted (mounting bracket included), set up on a table or on a single stand.

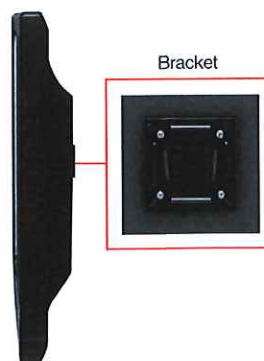


# EDC-2600

## Digital Chart LCD

### Variable Working Distance With Automatic Chart Calibration

The chart size automatically calibrates to the working distance, eliminating the need for mirror systems that compromise image quality, as well as the need for a technician. The working distance may be between 6' to 21' (2m - 6.5 m)

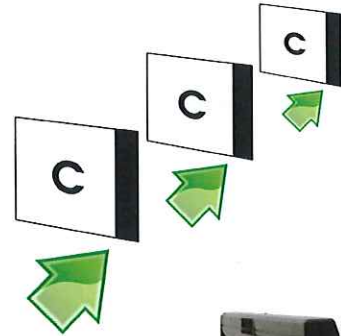


### Fast And Silent Picture Loading

Since the EDC-2600 software is based on a SD card, picture loading is extremely fast.



### Independent Software System

Upgrading your EDC-2600 is simple since the chart software is located in an SD card. When it is time for an upgrade, simply replace the SD card for an easy and fast process.

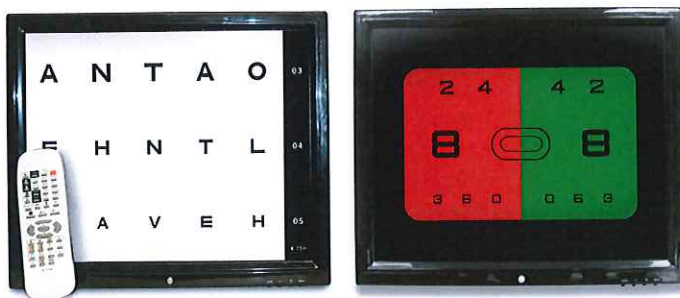

Snellen Acuity



ELI  
**EZER**

### Auto Saving Mode

The EDC-2600 features an auto saving mode that turns the monitor off after 10 minutes of inactivity.





# EDC-2600

Digital Chart LCD

## Specifications

<b>Panel</b>	19" (4:3)
<b>Pixel</b>	pitch 0.285mm
<b>Max brightness</b>	300cd/m2
<b>Max contrast ratio</b>	800:1
<b>Max resolution</b>	1280x1024
<b>Input signal</b>	RGB
<b>Signal joint</b>	15Pin D Pattern Plug
<b>Energy consumption</b>	35W
<b>Net weight</b>	5Kg
<b>Compatibility</b>	PC/MAC
<b>Dimension</b>	60(D) x 425(W) x 440(H)mm

